

Hox-A7 rabbit pAb**Cat#: orb765433 (Manual)**


For research use only. Not intended for diagnostic use.

Product Name	Hox-A7 rabbit pAb
Host species	Rabbit
Applications	WB;ELISA
Species Cross-Reactivity	Human;Mouse;Rat
Recommended dilutions	Western Blot: 1/500 - 1/2000. ELISA: 1/40000. Not yet tested in other applications.
Immunogen	The antiserum was produced against synthesized peptide derived from human HOXA7. AA range:91-140
Specificity	Hox-A7 Polyclonal Antibody detects endogenous levels of Hox-A7 protein.
Formulation	Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide..
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Protein Name	Homeobox protein Hox-A7
Gene Name	HOXA7
Cellular localization	Nucleus.
Purification	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.
Clonality	Polyclonal


Concentration	1 mg/ml
Observed band	25kD
Human Gene ID	3204
Human Swiss-Prot Number	P31268
Alternative Names	HOXA7; HOX1A; Homeobox protein Hox-A7; Homeobox protein Hox 1.1; Homeobox protein Hox-1A

Background

In vertebrates, the genes encoding the class of transcription factors called homeobox genes are found in clusters named A, B, C, and D on four separate chromosomes. Expression of these proteins is spatially and temporally regulated during embryonic development. This gene is part of the A cluster on chromosome 7 and encodes a DNA-binding transcription factor which may regulate gene expression, morphogenesis, and differentiation. For example, the encoded protein represses the transcription of differentiation-specific genes during keratinocyte proliferation, but this repression is then overcome by differentiation signals. This gene is highly similar to the antennapedia (Antp) gene of *Drosophila*. [provided by RefSeq, Jul 2008],

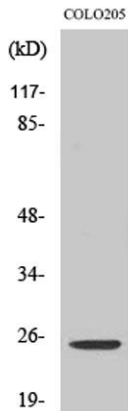

Western Blot analysis of various cells using Hox-A7 Polyclonal Antibody cells nucleus extracted by Minute TM Cytoplasmic and Nuclear Fractionation kit (SC-003, Invent biotech, MN, USA).

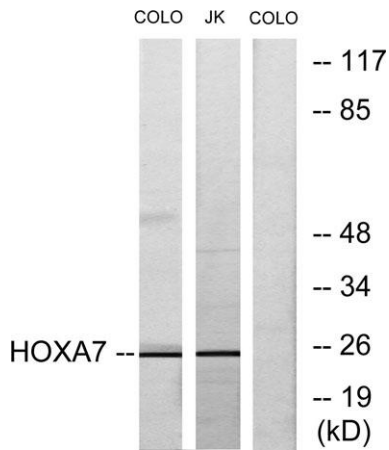

Western blot analysis of lysates from COLO205 and Jurkat cells, using HOXA7 Antibody. The lane on the right is blocked with the synthesized peptide.