

PAK4/5/6 rabbit pAb**Cat#: orb767333 (Manual)**


For research use only. Not intended for diagnostic use.

Product Name	PAK4/5/6 rabbit pAb
Host species	Rabbit
Applications	WB;ELISA
Species Cross-Reactivity	Human;Mouse;Rat
Recommended dilutions	Western Blot: 1/500 - 1/2000. ELISA: 1/20000. Not yet tested in other applications.
Immunogen	The antiserum was produced against synthesized peptide derived from human PAK4/5/6. AA range:441-490
Specificity	PAK4/5/6 Polyclonal Antibody detects endogenous levels of PAK4/5/6 protein.
Formulation	Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide..
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Protein Name	Serine/threonine-protein kinase PAK 4
Gene Name	PAK4 PAK5 PAK6
Cellular localization	Cytoplasm . Seems to shuttle between cytoplasmic compartments depending on the activating effector. For example, can be found on the cell periphery after activation of growth-factor or integrin-mediated signaling pathways. .
Purification	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.


Clonality	Polyclonal
Concentration	1 mg/ml
Observed band	80kD
Human Gene ID	10298
Human Swiss-Prot Number	O96013;Q9P286;Q9NQU5
Alternative Names	PAK4; KIAA1142; Serine/threonine-protein kinase PAK 4; p21-activated kinase 4; PAK-4

Background

PAK proteins, a family of serine/threonine p21-activating kinases, include PAK1, PAK2, PAK3 and PAK4. PAK proteins are critical effectors that link Rho GTPases to cytoskeleton reorganization and nuclear signaling. They serve as targets for the small GTP binding proteins Cdc42 and Rac and have been implicated in a wide range of biological activities. PAK4 interacts specifically with the GTP-bound form of Cdc42Hs and weakly activates the JNK family of MAP kinases. PAK4 is a mediator of filopodia formation and may play a role in the reorganization of the actin cytoskeleton. Multiple alternatively spliced transcript variants encoding distinct isoforms have been found for this gene. [provided by RefSeq, Jul 2008],

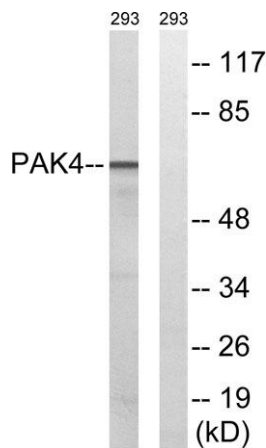

Western blot analysis of lysates from K562 cells, treated with PMA 125ng/ml 30', using PAK4/5/6 Antibody. The lane on the right is blocked with the synthesized peptide.

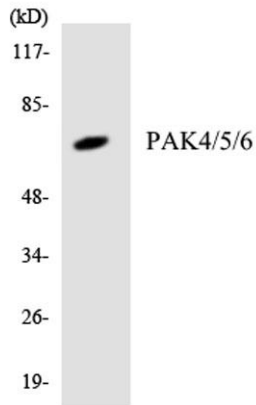

Western blot analysis of the lysates from Jurkat cells using PAK4/5/6 antibody.