

USP19 rabbit pAb**Cat#: orb767439 (Manual)**


For research use only. Not intended for diagnostic use.

| | |
|---------------------------------|--|
| Product Name | USP19 rabbit pAb |
| Host species | Rabbit |
| Applications | WB;IHC;IF;ELISA |
| Species Cross-Reactivity | Human;Mouse;Rat |
| Recommended dilutions | Western Blot: 1/500 - 1/2000. Immunohistochemistry: 1/100 - 1/300. ELISA: 1/20000. Not yet tested in other applications. |
| Immunogen | The antiserum was produced against synthesized peptide derived from human USP19. AA range:391-440 |
| Specificity | USP19 Polyclonal Antibody detects endogenous levels of USP19 protein. |
| Formulation | Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide.. |
| Storage | Store at -20°C. Avoid repeated freeze-thaw cycles. |
| Protein Name | Ubiquitin carboxyl-terminal hydrolase 19 |
| Gene Name | USP19 |
| Cellular localization | Endoplasmic reticulum membrane ; Single-pass membrane protein . |
| Purification | The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen. |
| Clonality | Polyclonal |


| | |
|--------------------------------|---|
| Concentration | 1 mg/ml |
| Observed band | 130kD |
| Human Gene ID | 10869 |
| Human Swiss-Prot Number | O94966 |
| Alternative Names | USP19; KIAA0891; ZMYND9; Ubiquitin carboxyl-terminal hydrolase 19; Deubiquitinating enzyme 19; Ubiquitin thioesterase 19; Ubiquitin-specific-processing protease 19; Zinc finger MYND domain-containing protein 9 |

Background


catalytic activity:Ubiquitin C-terminal thioester + H(2)O = ubiquitin + a thiol.,function:Involved in the ubiquitin-dependent proteolytic pathway in conjunction with the 26S proteasome.,PTM:Phosphorylated upon DNA damage, probably by ATM or ATR.,sequence caution:Translation N-terminally shortened.,similarity:Belongs to the peptidase C19 family.,similarity:Contains 1 MYND-type zinc finger.,similarity:Contains 2 CS domains.,

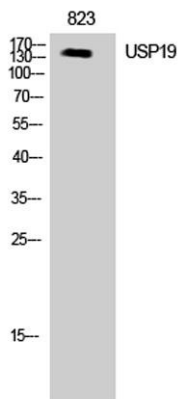

Western Blot analysis of NIH-3T3 cells using USP19 Polyclonal Antibody diluted at 1:2000. Secondary antibody(catalog#:RS0002) was diluted at 1:20000

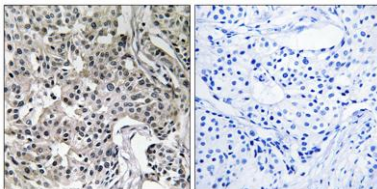

Immunohistochemistry analysis of paraffin-embedded human breast carcinoma tissue, using USP19 Antibody. The picture on the right is blocked with the synthesized peptide.

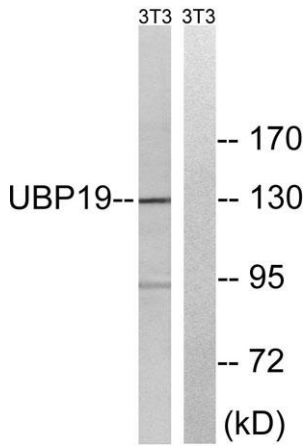

Western blot analysis of lysates from NIH/3T3 cells, using USP19 Antibody. The lane on the right is blocked with the synthesized peptide.