

**Versican rabbit pAb****Cat#: orb767741 (Manual)**


For research use only. Not intended for diagnostic use.

<b>Product Name</b>	Versican rabbit pAb
<b>Host species</b>	Rabbit
<b>Applications</b>	WB;IHC;IF;ELISA
<b>Species Cross-Reactivity</b>	Human;Rat;Mouse;
<b>Recommended dilutions</b>	Immunohistochemistry: 1/100 - 1/300. ELISA: 1/10000. Not yet tested in other applications.
<b>Immunogen</b>	The antiserum was produced against synthesized peptide derived from human VCAN. AA range:532-581
<b>Specificity</b>	Versican Polyclonal Antibody detects endogenous levels of Versican protein.
<b>Formulation</b>	Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide..
<b>Storage</b>	Store at -20°C. Avoid repeated freeze-thaw cycles.
<b>Protein Name</b>	Versican core protein
<b>Gene Name</b>	VCAN
<b>Cellular localization</b>	Secreted, extracellular space, extracellular matrix . Cell projection, cilium, photoreceptor outer segment . Secreted, extracellular space, extracellular matrix, interphotoreceptor matrix .
<b>Purification</b>	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.


<b>Clonality</b>	Polyclonal
<b>Concentration</b>	1 mg/ml
<b>Observed band</b>	
<b>Human Gene ID</b>	1462
<b>Human Swiss-Prot Number</b>	P13611
<b>Alternative Names</b>	VCAN; CSPG2; Versican core protein; Chondroitin sulfate proteoglycan core protein 2; Chondroitin sulfate proteoglycan 2; Glial hyaluronate-binding protein; GHAP; Large fibroblast proteoglycan; PG-M

## Background

This gene is a member of the aggrecan/versican proteoglycan family. The protein encoded is a large chondroitin sulfate proteoglycan and is a major component of the extracellular matrix. This protein is involved in cell adhesion, proliferation, proliferation, migration and angiogenesis and plays a central role in tissue morphogenesis and maintenance. Mutations in this gene are the cause of Wagner syndrome type 1. Multiple transcript variants encoding different isoforms have been found for this gene. [provided by RefSeq, Aug 2009],

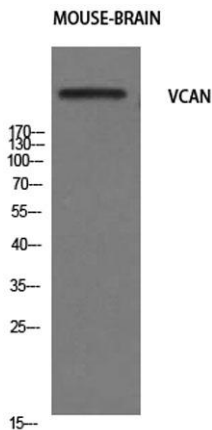

**Western Blot analysis of mouse-brain cells using Versican Polyclonal Antibody. Secondary antibody(catalog#:RS0002) was diluted at 1:20000**

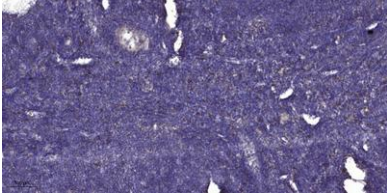

**Immunohistochemical analysis of paraffin-embedded human tonsil. 1, Antibody was diluted at 1:200(4° overnight). 2, Tris-EDTA,pH9.0 was used for antigen retrieval. 3,Secondary antibody was diluted at 1:200(room temperature, 30min).**