

PKA II β reg (phospho Ser113) rabbit pAb**Cat#: orb769634 (Manual)**


For research use only. Not intended for diagnostic use.

Product Name	PKA II β reg (phospho Ser113) rabbit pAb
Host species	Rabbit
Applications	WB;IHC;IF;ELISA
Species Cross-Reactivity	Human;Mouse;Rat;Monkey
Recommended dilutions	Western Blot: 1/500 - 1/2000. Immunohistochemistry: 1/100 - 1/300. ELISA: 1/10000. Not yet tested in other applications.
Immunogen	The antiserum was produced against synthesized peptide derived from human PKA-R2 beta around the phosphorylation site of Ser113. AA range:79-128
Specificity	Phospho-PKA II β reg (S113) Polyclonal Antibody detects endogenous levels of PKA II β reg protein only when phosphorylated at S113.
Formulation	Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide..
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Protein Name	cAMP-dependent protein kinase type II-beta regulatory subunit
Gene Name	PRKAR2B
Cellular localization	Cytoplasm . Cell membrane . Colocalizes with PJA2 in the cytoplasm and at the cell membrane.
Purification	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.
Clonality	Polyclonal


Concentration	1 mg/ml
Observed band	46kD
Human Gene ID	5577
Human Swiss-Prot Number	P31323
Alternative Names	PRKAR2B; cAMP-dependent protein kinase type II-beta regulatory subunit

Background

cAMP is a signaling molecule important for a variety of cellular functions. cAMP exerts its effects by activating the cAMP-dependent protein kinase, which transduces the signal through phosphorylation of different target proteins. The inactive kinase holoenzyme is a tetramer composed of two regulatory and two catalytic subunits. cAMP causes the dissociation of the inactive holoenzyme into a dimer of regulatory subunits bound to four cAMP and two free monomeric catalytic subunits. Four different regulatory subunits and three catalytic subunits have been identified in humans. The protein encoded by this gene is one of the regulatory subunits. This subunit can be phosphorylated by the activated catalytic subunit. This subunit has been shown to interact with and suppress the transcriptional activity of the cAMP responsive element binding protein 1 (CREB1) in activ

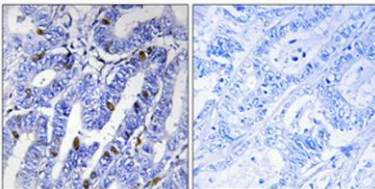

Immunohistochemical analysis of paraffin-embedded Human colon cancer. Antibody was diluted at 1:100(4° overnight). High-pressure and temperature Tris-EDTA,pH8.0 was used for antigen retrieval. Negative contrl (right) obtained from antibody was pre-absorb

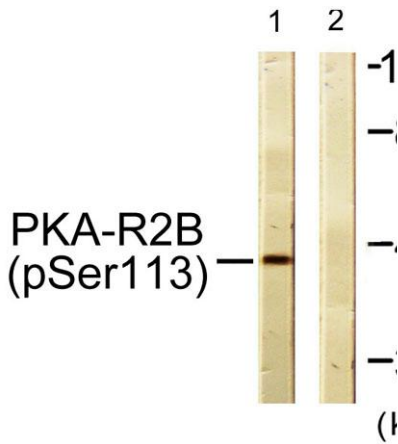

Western blot analysis of lysates from COS7 cells treated with PMA 125ng/ml 30', using PKA-R2 beta (Phospho-Ser113) Antibody. The lane on the right is blocked with the phospho peptide.