

RUNX1 (phospho Ser249) rabbit pAb**Cat#: orb770708 (Manual)**

For research use only. Not intended for diagnostic use.


Product Name	RUNX1 (phospho Ser249) rabbit pAb
Host species	Rabbit
Applications	WB;ELISA
Species Cross-Reactivity	Human;Mouse;Rat
Recommended dilutions	Western Blot: 1/500 - 1/2000. ELISA: 1/5000. Not yet tested in other applications.
Immunogen	The antiserum was produced against synthesized peptide derived from human AML1 around the phosphorylation site of Ser276. AA range:242-291
Specificity	Phospho-RUNX1 (S249) Polyclonal Antibody detects endogenous levels of RUNX1 protein only when phosphorylated at S249.
Formulation	Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide..
Storage	Store at -20°C. Avoid repeated freeze-thaw cycles.
Protein Name	Runt-related transcription factor 1
Gene Name	RUNX1
Cellular localization	Nucleus.
Purification	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.
Clonality	Polyclonal

Explore. Bioreagents.


Concentration	1 mg/ml
Observed band	49kD
Human Gene ID	861
Human Swiss-Prot Number	Q01196
Alternative Names	RUNX1; AML1; CBFA2; Runt-related transcription factor 1; Acute myeloid leukemia 1 protein; Core-binding factor subunit alpha-2; CBF-alpha-2; Oncogene AML-1; Polyomavirus enhancer-binding protein 2 alpha B subunit; PEA2-alpha B; PEBP2-alpha

Background

Core binding factor (CBF) is a heterodimeric transcription factor that binds to the core element of many enhancers and promoters. The protein encoded by this gene represents the alpha subunit of CBF and is thought to be involved in the development of normal hematopoiesis. Chromosomal translocations involving this gene are well-documented and have been associated with several types of leukemia. Three transcript variants encoding different isoforms have been found for this gene. [provided by RefSeq, Jul 2008],

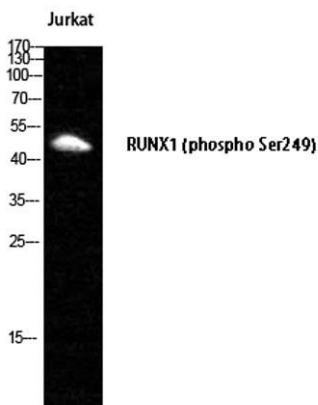

Western Blot analysis of Jurkat cells using Phospho-RUNX1 (S249) Polyclonal Antibody cells nucleus extracted by Minute TM Cytoplasmic and Nuclear Fractionation kit (SC-003, Inventbiotech, MN, USA).

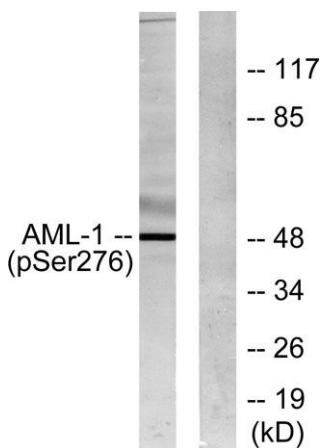

Western blot analysis of lysates from Jurkat cells, using AML1 (Phospho-Ser276) Antibody. The lane on the right is blocked with the phospho peptide.



Explore. Bioreagents.

www.biorbyt.com