

---

## Product Datasheet

### ABO Antibody (orb1252539)

**Description**

ABO Antibody

**Species/Host**

Mouse

**Reactivity**

Human

**Conjugation**

Unconjugated

**Tested Applications**

FA, IF, IHC-P, WB

**Immunogen**

A mixture of erythrocytes of group B and glycoprotein fraction isolated from saliva of secretors with blood group B was used as the immunogen for the ABO antibody.

**Target**

ABO

**Preservatives**

PBS with 0.1 mg/ml rAlbumin and 0.05% sodium azide

**Form/Appearance**

Liquid

**Concentration**

0.2 mg/mL

**Storage**

Aliquot and Store at 2-8°C. Avoid freeze-thaw cycles.

**Note**

For research use only

**Application notes**

Functional Activity-Erythrocyte Agglutination  
 Immunofluorescence: 0.5-1 ug/ml  
 Immunohistochemistry (FFPE): 0.5-1 ug/ml for 30 min at RT  
 (1) Prediluted format : incubate for 30 min at RT  
 (2) Optimal dilution of the ABO antibody should be determined by the researcher.  
 1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes  
 2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

**Isotype**


IgM, kappa

**Clonality**


Monoclonal

**Uniprot ID**
[P16442](#)
**Dilution Range**

Functional Activity-Erythrocyte Agglutination  
 Immunofluorescence: 0.5-1 ug/ml  
 Immunohistochemistry (FFPE): 0.5-1 ug/ml for 30 min at RT  
 (1) Prediluted format : incubate for 30 min at RT  
 (2) Optimal dilution of the ABO antibody should be determined by the researcher.  
 1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes  
 2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if

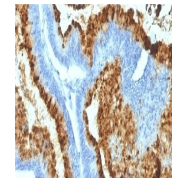

Formalin-fixed, paraffin-embedded human ...

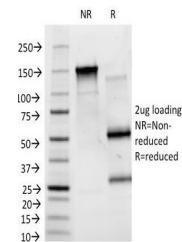

SDS-PAGE Analysis of Purified, BSA-Free ...