

## **Product Datasheet**

## Cattle PSPB(Pregnancy Specific Protein B) ELISA Kit (orb1817306)

**Description** 

This assay employs the competitive inhibition enzyme immunoassay technique. The microtiter plate provided in this kit has been pre-coated with Cattle PSPB protein. Standards or samples are added to the appropriate microtiter plate wells then with a biotin-conjugated antibody specific to Cattle PSPB. Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of  $450 \, \mathrm{nm} \pm 10 \, \mathrm{nm}$ . The concentration of Cattle PSPB in the samples is then determined by comparing the OD of the samples to the standard curve.

**Reactivity** Bovine

**Range** 0.78-50 ng/mL

**Concentration** 50 ng/mL

**Note** For research use only

**Application notes** 

standard: 50 ng/mL. Test principle: This assay employs the competitive inhibition enzyme immunoassay technique. The microtiter plate provided in this kit has been pre-coated with Cattle PSPB protein. Standards or samples are added to the appropriate microtiter plate wells then with a biotin-conjugated antibody specific to Cattle PSPB. Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of 450nm  $\pm$  10nm. The concentration of Cattle PSPB in the samples is then determined by comparing the OD of the samples to the standard curve


**Sample Types** serum, plasma and other biological fluids

Sensitivity 0.29 ng/mL

**Expiration Date** Please enquire.





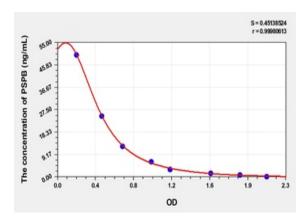

## **Biorbyt LLC.**

Email: <u>info@biorbyt.com</u>, <u>support@biorbyt.com</u>
Phone: <u>+1 (415) 906-5211</u> | Fax: <u>+1 (415) 651-8558</u>