

Product Datasheet

Goat Anti-Mouse/Rabbit IgG (H&L)-HRP Polymer (Ready to use) (orb1878754)


Description	Ready to use Goat Polyclonal Secondary Antibody to Mouse/Rabbit IgG (H&L) HRP polymer labeled
Species/Host	Goat
Reactivity	Mouse, Rabbit
Conjugation	HRP
Tested Applications	ELISA, IH, WB
Immunogen	Mouse/Rabbit IgG
Form/Appearance	Liquid in PBS, pH 7.2, containing 1% rAlbumin, 0.05% Proclin
Storage	Shipped and store at 4°C for one year.
Note	For research use only
Clonality	Polyclonal
Source	Goat
Dilution Range	Ready to use
Expiration Date	12 months from date of receipt.

Biorbyt Ltd.


7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Biorbyt LLC.

68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+1(415)9065211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

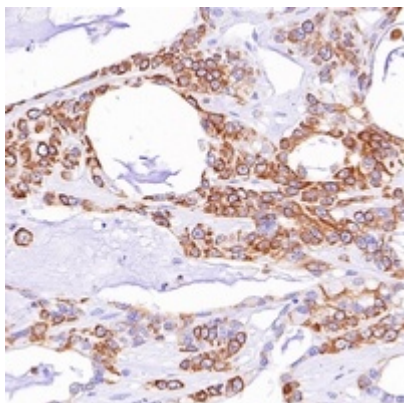

Immunohistochemical analysis staining in human liver carcinoma formalin fixed paraffin-embedded tissue section. The section was pre-treated using pressure cooker heat antigen retrieval with sodium citrate buffer (0.01 M, pH = 6) for 3 minutes. The section was detected using mouse primary antibody, and Goat Anti-Mouse/Rabbit IgG (H&L)-HRP polymer (Ready to use). The section was then counterstained with haematoxylin and mounted with Neutral Gum.

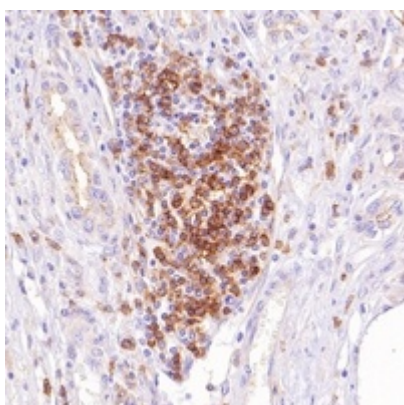

Immunohistochemical analysis staining in human liver carcinoma formalin fixed paraffin-embedded tissue section. The section was pre-treated using pressure cooker heat antigen retrieval with sodium citrate buffer (0.01 M, pH = 6) for 3 minutes. The section was detected using rabbit primary antibody, and Goat Anti-Mouse/Rabbit IgG (H&L)-HRP polymer (Ready to use). The section was then counterstained with haematoxylin and mounted with Neutral Gum.

Biorbyt Ltd.

7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Biorbyt LLC.

68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+1(415)9065211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)