

Product Datasheet

Guinea Pig IgG (H&L) Antibody Texas Red Conjugated Pre-Adsorbed (orb347056)

Description Guinea Pig IgG (H&L) antibody (Texas Red)

Species/Host Goat

Reactivity Guinea pig

Conjugation Texas Red

Tested Applications FC, FLISA, IF

Immunogen Guinea Pig IgG whole molecule

Preservatives 0.01% (w/v) Sodium Azide

Form/Appearance Lyophilized

Concentration 1.0 mg/mL

Storage Store vial at 4° C prior to restoration. For extended storage aliquot contents and

freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to

immediate use.

Note For research use only

Application notes This product is designed for immunofluorescence microscopy, fluorescence

based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various

commercial platforms.

Isotype lgG

Clonality Polyclonal





Purity


This product was prepared from monospecific antiserum by immunoaffinity chromatography using Guinea Pig IgG coupled to agarose beads followed by solid phase adsorption(s) to remove any unwanted reactivities. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Goat Serum, Guinea Pig IgG and Guinea Pig Serum. No reaction was observed against Bovine, Chicken, Goat, Hamster, Horse, Human, Mouse, Rabbit, Rat and Sheep Serum Proteins.

Dilution Range


FLISA: 1:10,000 - 1:50,000, FC: 1:500 - 1:2,500, IF: 1:1,000 - 1:5,000

Expiration Date

12 months from date of receipt.

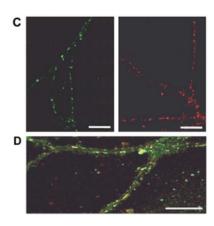


(C) Surface distribution of anti-GluR1-4R (green) and anti-GluR1-4GP (red) immunoreactivity on living hippocampal neurons in culture (14 DIV). (D) Colocalization (yellow) of antiGluR1-4R (green) and anti-GluR1-4GP (red) immunoreactivity on the surface of hippocampal neurons. Neuronal cultures were first stained with the anti-GluR1-4R (green) antibody, followed by restaining with the anti-GluR1-4GP (red) antibody under non-permeabilised conditions (as described at Fig. 1A, without the LTP induction). Scale bars 20 μm .

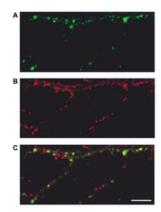

Changes in AMPAR expression on the surface of living hippocampal neurons following high K + treatment. Dual labelling of hippocampal dendrites showing AMPAR clusters labelled before (A; green and yellow) and following (B; red) the high K + treatment. (C) Colocalization (yellow) of anti-GluR1-4R (green) and anti-GluR1-4GP (red) immunoreactivity on the surface of hippocampal neurons. Note the marked increase in new AMPAR clusters (C; red only). Scale bar 5 µm.





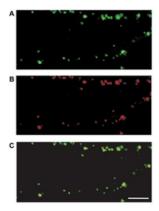

Dual labelling of hippocampal dendrites before (A; green) and following (B; red) the high K+ treatment in the presence of an NMDAR antagonist (5 μ m L-689, 560). Colocalization (yellow) of anti-GluR1-4R (green) and anti-GluR1-4GP (red) immunoreactivity on the surface of hippocampal neurons. Note the absence of new (red only) AMPAR clusters on panel C. Scale bar 5 μ m.

Phone: +1 (415) 906-5211 | Fax: +1 (415) 651-8558