

Product Datasheet

Pan Cytokeratin Antibody Cocktail (Acidic + Basic) (orb639459)

Catalog Number	orb639459
Category	Antibodies
Description	Twenty human keratins are resolved with two-dimensional gel electrophoresis into acidic (pI 6.0) subfamilies. This pan keratin antibody cocktail recognizes acidic (Type I or LMW) and basic (Type II or HMW) cytokeratins, which include CK1, CK3-6, CK8, CK10, CK14-16, and CK19. Many studies have shown the usefulness of keratin markers in cancer research and tumor diagnosis. The AE1 + AE3 antibody cocktail is a broad spectrum pan cytokeratin antibody cocktail which differentiates epithelial tumors from non-epithelial tumors e.g. squamous vs. adenocarcinoma of the lung, liver carcinoma, breast cancer, and esophageal cancer. It has been used to characterize the source of various neoplasms and to study the distribution of keratin containing cells in epithelia during normal development and during the development of epithelial neoplasms. It stains cytokeratin present in normal and abnormal human tissues and has shown high sensitivity in the recognition of epithelial cells and carcinomas.
Clonality	Monoclonal
Species/Host	Mouse
Isotype	Mouse IgG1, kappa
Conjugation	Unconjugated
Reactivity	Human, Mouse, Rat
Immunogen	Human epidermal keratin was used as the immunogen for this pan Cytokeratin antibody.
Tested applications	FACS, IF, IHC, WB

Biorbyt Ltd.

7 Signet Court, Swann Road
Cambridge
CB5 8LA
United Kingdom


Email: info@biorbyt.com, support@biorbyt.com
Phone: +44 (0)1223 859353 | Fax: +1 (415) 651-8558

Biorbyt LLC

68 TW Alexander Drive
Research Triangle Park
Durham
NC 27713
United States

Email: info@biorbyt.com, support@biorbyt.com
Phone: +1 (415) 906-5211 | Fax: +1 (415) 651-8558

Dilution range	FACS: 0.5-1ug/million cells,Immunofluorescence: 1-2ug/ml,Western blot: 0.5-1ug/ml for 2 hours at RT,Immunohistochemistry (FFPE): 0.5-1ug/ml for 30 min at RT (1),Prediluted IHC only format: incubate for 30 min at RT (2)
Application notes	The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the pan Cytokeratin antibody AE1 + AE3 to be titrated up or down for optimal performance.1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes.2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.
Antibody Type	Primary Antibody
Clone Number	AE1 + AE3
Formula	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide
Storage	Maintain refrigerated at 2-8°C for up to 2 weeks. For long term storage store at -20°C in small aliquots to prevent freeze-thaw cycles.
Note	For research use only

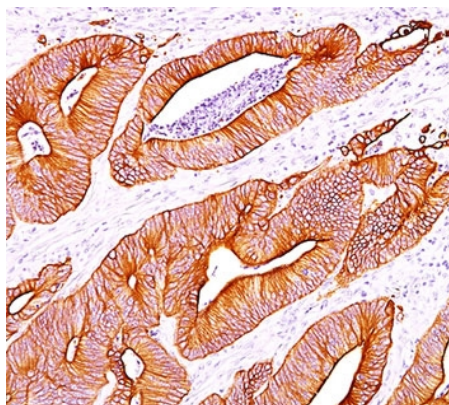

IHC staining of colon carcinoma with pan Cytokeratin antibody cocktail AE1 + AE3.

Biorbyt Ltd.


7 Signet Court, Swann Road
Cambridge
CB5 8LA
United Kingdom

Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\)1223 859353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)651-8558)


Biorbyt LLC

68 TW Alexander Drive
Research Triangle Park
Durham
NC 27713
United States

Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+1(415)906-5211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)651-8558)

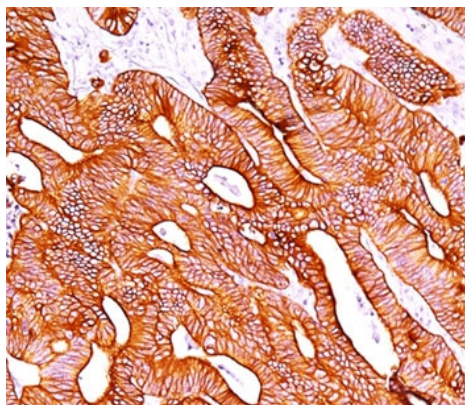

IHC staining of colon carcinoma with pan Cytokeratin antibody cocktail AE1 + AE3.

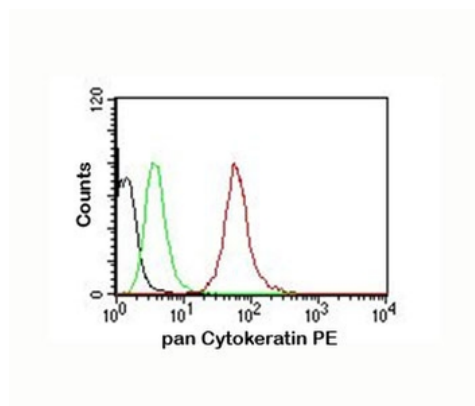

FACS testing of MCF-7 cells: Black = cells alone; Green = isotype control; Red = pan Cytokeratin antibody PE conjugate

Biorbyt Ltd.

7 Signet Court, Swann Road
Cambridge
CB5 8LA
United Kingdom

Email: info@biorbyt.com, support@biorbyt.com

Phone: [+44 \(0\)1223 859353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Biorbyt LLC

68 TW Alexander Drive
Research Triangle Park
Durham
NC 27713
United States

Email: info@biorbyt.com, support@biorbyt.com

Phone: [+1 \(415\) 906-5211](tel:+1(415)9065211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)