

Product Datasheet

CD79a Antibody (orb749344)

Description	A disulphide-linked heterodimer, consisting of mb-1 (or CD79a) and B29 (or CD79b) polypeptides, is non-covalently associated with membrane-bound immunoglobulins on B cells. This complex of mb-1 and B29 polypeptides and immunoglobulin constitute the B-cell Ag receptor. CD79a first appears at pre B-cell stage, early in maturation, and persists until the plasma cell stage where it is found as an intracellular component. It is found in the majority of acute leukemias of B-cell precursors, lines and lymphomas, and in some myelomas. It is not present in myeloid or Tcell lines. CD79a antibody is generally used to complement CD20 mAb, especially for mature B-cell lymphomas after treatment with Rituximab (anti-CD20). This CD79a antibody will stain many of the same lymphomas as CD20 mAb, but also is more likely to stain B-lymphoblastic lymphoma/leukemia than is CD20 mAb. CD79a antibody also stains more cases of plasma cell myeloma and occasionally some types of endothelial cells as well.
Species/Host	Mouse
Reactivity	Human
Conjugation	Unconjugated
Tested Applications	FACS, IF, IHC-P, WB
Immunogen	A soluble form of the extracellular IgSF domain from the human protein was used as the immunogen for the CD79a antibody.
Preservatives	0.2 mg/ml in 1X PBS with 0.1 mg/ml rAlbumin (US sourced) and 0.05% sodium azide
Storage	Store the CD79a antibody at 2-8°C (with azide) or aliquot and store at -20°C or colder (without azide).
Note	For research use only


Biorbyt Ltd.

7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Biorbyt LLC.

68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+1(415)9065211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Application notes	The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the CD79a antibody to be titered up or down for optimal performance.1. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.
Formula	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide
Isotype	Mouse IgG1, kappa
Clonality	Monoclonal
Clone Number	JCB117
Purity	Protein G affinity chromatography
Hazard Information	This CD79a antibody is available for research use only.
Entrez	973
Dilution Range	Western blot: 1-2ug/ml,Flow cytometry: 1-2ug/million cells,Immunofluorescence: 1-2ug/ml,Immunohistochemistry (FFPE): 1-2ug/ml for 30 min at RT
Expiration Date	12 months from date of receipt.

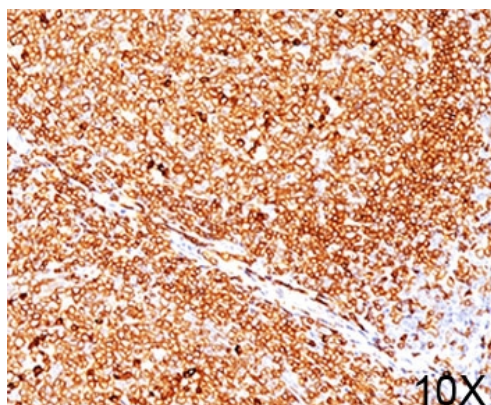


IHC testing of FFPE human tonsil stained with CD79a antibody (clone JCB117). Required HIER: boil tissue sections in pH9 10mM Tris with 1mM EDTA for 10-20 min followed by cooling at RT for 20 minutes.

Biorbyt Ltd.


7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Biorbyt LLC.


68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+1(415)9065211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

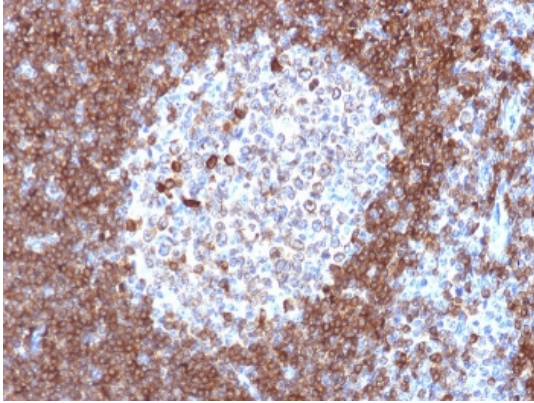

IHC testing of FFPE human tonsil stained with CD79a antibody (clone JCB117). Required HIER: boil tissue sections in pH9 10mM Tris with 1mM EDTA for 10-20 min followed by cooling at RT for 20 minutes.

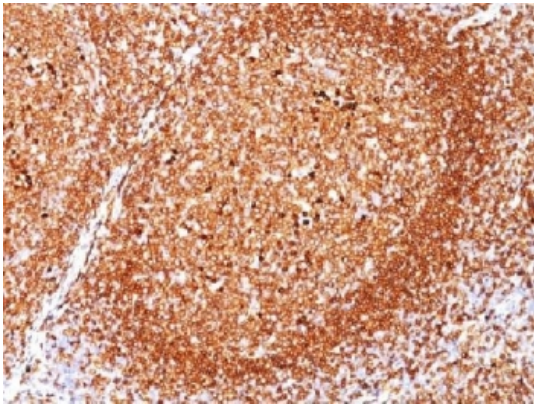

IHC testing of FFPE human tonsil stained with CD79a antibody (clone JCB117). Required HIER: boil tissue sections in pH9 10mM Tris with 1mM EDTA for 10-20 min followed by cooling at RT for 20 minutes.

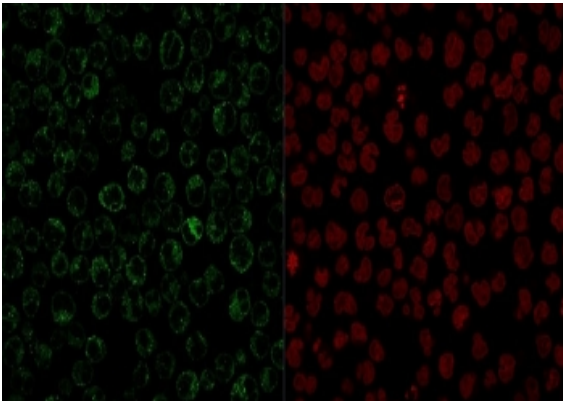


Immunofluorescent staining of PFA-fixed human Raji cells with CD79 antibody (clone JCB117, green) and Reddot nuclear stain (red).

Biorbyt Ltd.


7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Biorbyt LLC.


68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+1(415)9065211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

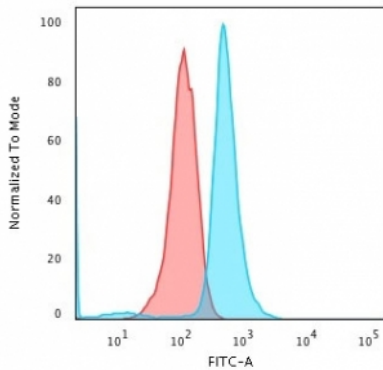

Flow cytometry testing of human Raji cells with CD79a antibody (clone JCB117); Red = isotype control, Blue = CD79a antibody.

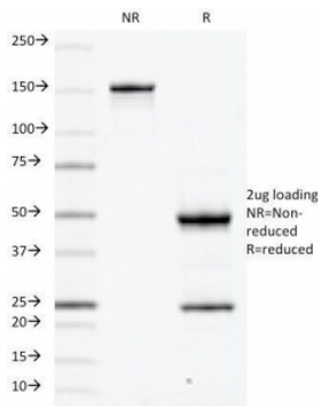

SDS-PAGE analysis of purified, BSA-free CD79a antibody (clone JCB117) as confirmation of integrity and purity.

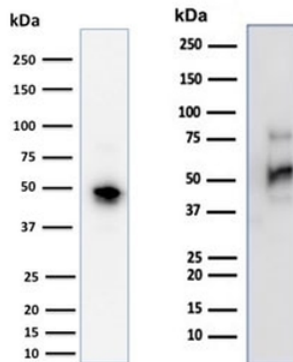

Western blot testing of two different lots of human Raji cell lysate with CD79a antibody (clone JCB117). Expected molecular weight: 25-47 kDa depending on glycosylation level.

Biorbyt Ltd.

7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Biorbyt LLC.

68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+1(415)9065211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)